



Esco Powdermax™ 1 Powder Weighing Balance Enclosure
Model PW1-3A_

Powdermax™

Powder Weighing Balance Enclosure

Introduction

Laboratory operations involving weighing, mixing and dispensing hazardous powders can now be performed safely and economically with the Esco Powdermax™ 1 Powder Weighing Balance Enclosure. The enclosure's unique design maintains balance stability to 4 decimal places without compromising containment.

Designed specifically for powder containment applications, the Powdermax Enclosure employs carefully directed inward airflow at the face of the enclosure to protect the operator from exposure to hazardous powders or potent pharmaceutical compounds. The aerodynamically designed enclosure reduces airflow turbulence for maximum containment while minimizing the effect that air drafts may have on sensitive weighing balances. The enclosure design exceeds the containment requirements of ANSI/ASHRAE 110 as proven by independent type testing. All exhaust air is HEPA-filtered, removing virtually all particulate matter, before recirculation to the laboratory.

Designed for Enhanced Usability and Efficiency

The Powdermax™ is designed specifically to provide the operator with a high level of usability, comfort and visibility.

- A 13-degree sloped front allows easy access to the work zone, eliminating operator fatigue and increasing productivity.
- Transparent 6 mm (0.2") frameless acrylic front and sides enhance visibility and operator comfort during longer periods of operation. The self-supporting pre-tensioned hinges on the front window provide easy access during loading and start-up.
- Ergonomically designed apertures for hands simultaneously provide maximum movement within the work zone while protecting the operator. The curved front edge minimizes airflow turbulence and improves user comfort.

- Electronic ballast for the fluorescent lighting provides zero-flicker with increased energy efficiency, reliability and service life with a lower heat output.
- Two electrical pass throughs on the back wall of the unit provide convenient access to power sources for equipment.
- The spacious interior dimensions will easily accommodate large analytical balances.

Enhanced Filtration System

The enclosure's unique filtration system ensures the highest level of protection for the operator.


- The Powdermax™ provides enhanced safety to laboratory personnel by drawing air across the front opening in the hood at a controlled rate, preventing the escape of powders and particulate matters into the general laboratory environment.
- A unique filter mounting system eliminates bypass leaks inherent in conventional clamping systems.
- The Powdermax™ is equipped with a generously sized HEPA filter system, rated at >99.99% efficiency for 0.3 micron sized particles.
- The aerodynamically designed system constructed of 1.2 mm (18 gauge) epoxy powder-coated electrogalvanized steel delivers maximum containment by ensuring airflow uniformity throughout the work chamber of the cabinet.
- Since HEPA-filtered air is recirculated to the laboratory, the enclosure does not have to be connected to an extraction/exhaust system, resulting in low operating costs and mobility.

User-Friendly Control System

The user-friendly Esco Sentinel™ Silver microprocessor based control system is fully configurable according to operator's requirements and comes equipped with a number of enhanced features to promote cabinet usability.

ESCO

WORLD CLASS. WORLDWIDE.

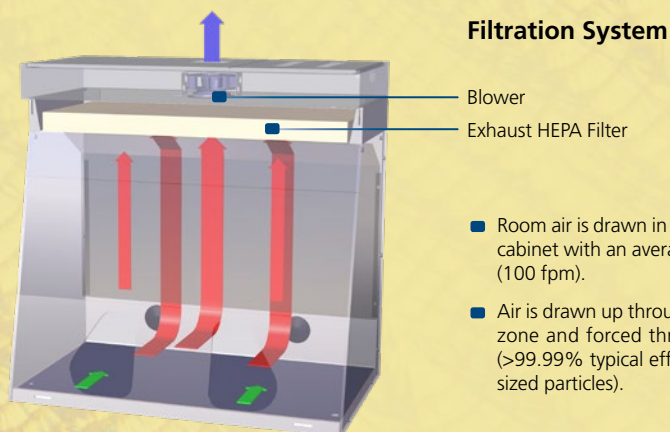

- HEPA-filtered air
- Unfiltered / potentially contaminated air
- Inflow air

- Room air is drawn in across the front of the cabinet with an average velocity of 0.5 m/s (100 fpm).
- Air is drawn up through the cabinet's work zone and forced through the HEPA filter (>99.99% typical efficiency for 0.3 micron sized particles).
- The state-of-the-art baffle system ensures airflow uniformity throughout the cabinet's main chamber.
- The HEPA filtered air then returns to the laboratory stripped of all airborne contaminants or is vented through the optional exhaust collar to exhaust ducting for enhanced safety.

- Accurate true airflow velocity sensing technology measures all critical cabinet airflow parameters allowing superior monitoring. Temperature compensated sensors ensure increased accuracy.
- Password-protected administration can be set to restrict access to the main menu.
- Solid state variable speed controllers offer superior control over conventional "step" controllers. The built in RFI and noise filters eliminate interference with adjacent instrumentation.
- A bright, easy-to-read, LCD display provides continuous monitoring of cabinet airflow.
- Audible and visual alarms ensure operator and environmental protection by alerting the user in the event of low airflow or unsafe sash positions.
- When the cabinet is switched on, its automatic warm-up cycle is activated, purging all contaminants from the work zone prior to use.
- An automatic post-purge cycle can be configured on shutdown, purging all

residual contaminants before deactivation.

Highest Quality Enclosure Construction

All Esco cabinets are constructed to the highest quality using the finest materials.

- The enclosure's head unit is constructed of electrogalvanized steel with abrasion-resistant powder-coated finish.
- External surfaces are coated with Esco ISOCIDE™ antimicrobial coating to protect against surface contamination and inhibit bacterial growth. Isocide eliminates 99.9% of surface bacteria within 24 hours of exposure.
- The enclosure is supplied with a chemical and abrasion resistant phenolic resin removable work surface. The solid work surface minimizes transmission of fan vibration to the user's weighing balance.

Blower Efficiency

The Powdermax™ incorporates an energy efficient external rotor motor design which reduces operating costs and has extremely low noise and vibration levels. The cabinet

blower system creates an inflow of air from the ambient environment into the cabinet with an average velocity of 0.5 m/s (100 fpm).

Designed and Built to Exceed Safety Criteria

All components used in Esco products meet or exceed all the applicable safety requirements and are UL listed/recognized. The Powdermax™ is compliant with International Standards. The cabinet is designed to meet the general safety requirements of EN 61010-1, EN 61326 Class B.

Options and Accessories

The Powdermax™ is available with a number of options and accessories to meet your needs. These include:

- Mobile cart with foldable trays
- Height-adjustable lab chair
- Ergonomic foot rest

Further customization specifications are available upon request.

Microbalance Stability Testing

Weighing Accuracy Test on Esco Powdermax™ Powder Weighing Balance Enclosure



The Esco PowderMax™ cabinet at nominal airflow can provide a stable work condition for sensitive powder weighing work with extremely low standard deviation in the range of 0.00001 to 0.000011 for weighing 20g, 5g, and even 0.5g objects.

Surrogate Powder Testing


Active Pharmaceutical Ingredients (API) Surrogate Containment Testing on Esco Powdermax™ Powder Weighing Balance Enclosure



The Esco PowderMax™ Powder Weighing Balance Enclosure, operating at 0.50 m/s (100fpm) average inflow, was demonstrated to contain surrogate API (naproxen sodium), at an OEL of 26.3 nanogram / m³, for a challenge of repetitive weighing of 20 x 1-gram quantities of API, independent of the operator's level of skill.

General Specifications, Powder Weighing Balance Enclosure			
Model	220-240 VAC, 50 Hz	PW-3A1 2160001	PW-4A1 2160010
	110-120 VAC, 60 Hz	PW-3A2 2160003	PW-4A2 2160006
	220-240 VAC, 60 Hz	PW-3A3 2160004	PW-4A3 2160007
Nominal Size		0.9 meters (3 ft)	1.2 meters (4 ft)
External Dimensions (W x D x H)		840 x 700 x 1125 mm (33.0" x 27.6" x 44.3")	1140 x 700 x 1125 mm (44.9" x 27.6" x 44.3")
Internal Work Area, Dimensions (W x D x H)		820 x 650 x 829 mm (32.3" x 25.6" x 32.6")	1120 x 650 x 848 mm (44.1" x 25.6" x 33.4")
Standard Filtration Elements	Main Filter	HEPA filter, typical efficiency of >99.99% at 0.3 microns, removes particles and aerosols	
Control System		Sentinel™ Silver Microprocessor Controller	
Inflow Air Velocity		Initial setpoint: average of 0.5 m/s (100 fpm)	
Sound Emission		<55 dBA at initial blower speed setting measured at typical operator work position	
Cabinet Construction	Main Body	1.2 mm (0.05") 18 gauge electrogalvanized steel with white oven-baked epoxy-polyester ISOCIDE™ antimicrobial powder coated finish	
	Back Wall		
	Front Window	6 mm (0.2") Acrylic	
	Side Walls		
	Work Top	Phenolic Resin	
Fluorescent Lamp Intensity at Zero Ambient		>560 lux (>52 foot-candles)	>550 lux (>51 foot-candles)
Sash Specification	Sloped Front	13°	
	Sash Type	Hinged	
Net Weight		110 kg (242.5 lbs)	160 kg (352.7 lbs)
Shipping Weight		120 kg (264.6 lbs)	170 kg (374.8 lbs)
Shipping Dimensions, Maximum (W x D x H)		1050 x 1100 x 720 mm (2314.9" x 2425.1" x 1587.33")	1300 x 1100 x 720 mm (2866.0" x 2425.1" x 1587.33")
Shipping Volume, Maximum		0.83 m³ (29.3 cu.ft.)	1.0 m³ (35.3 cu.ft.)


Powdermax™ Powder Weighing Balance Enclosure Engineering Drawing

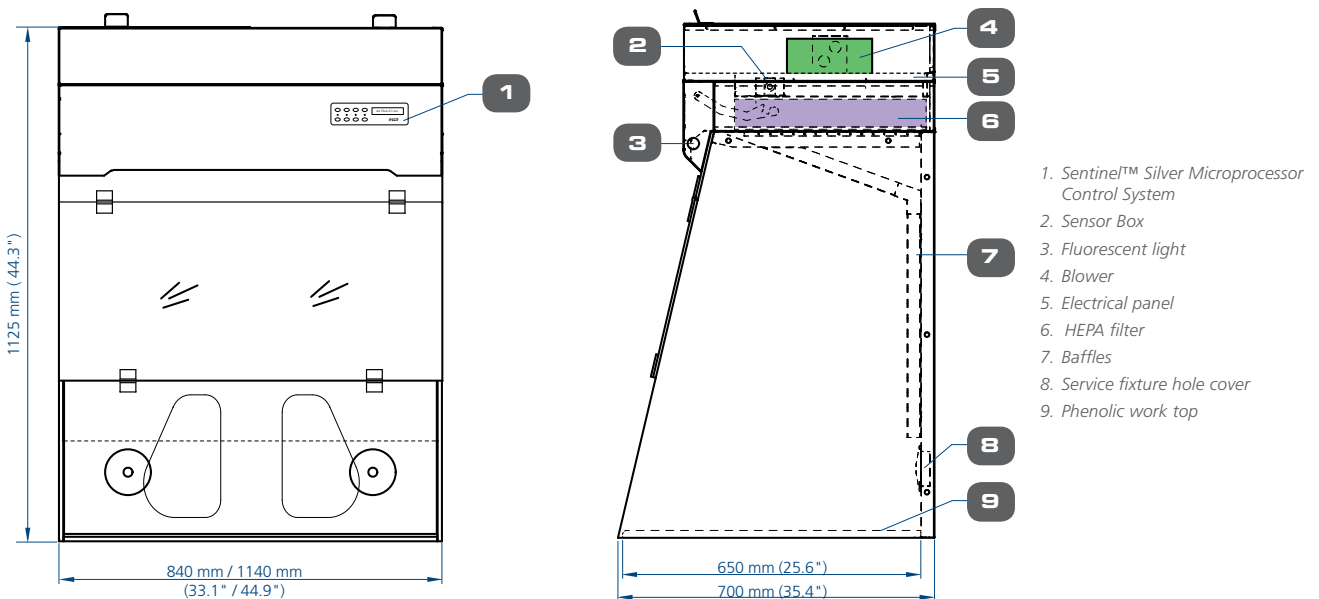

Standards Compliance	Containment	Electrical Safety
	ANSI/ASHRAE 110-1995, USA	UL-61010-1, USA CAN/CSA-C22-2 No.61010-1 EN-61010-1, Europe IEC61010-1, Worldwide

ESCO GLOBAL NETWORK

43 LOCATIONS IN 24 COUNTRIES ALL OVER THE WORLD

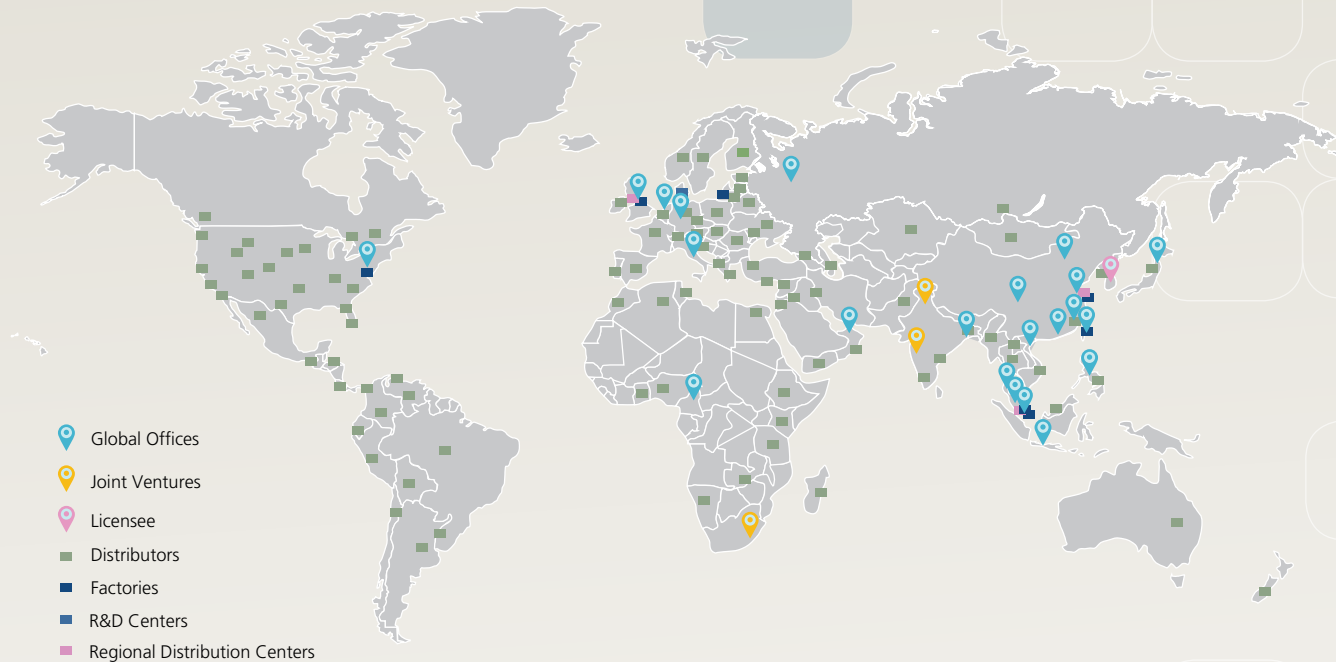

ART Equipment
 Biological Safety Cabinets
 CO₂ Incubators
 Compounding Pharmacy Equipment
 Containment / Pharma Products
 Ductless Fume Hoods
 Lab Animal Research Products
 Laboratory Centrifuges
 Laboratory Fume Hoods
 Laboratory Ovens and Incubators
 Laboratory Shakers
 Laminar Flow Clean Benches
 PCR Cabinets
 PCR Thermal Cyclers
 Powder Weighing Balance Enclosures
 Ultra-low Temperature Freezers

The Esco Group of Companies is a global life sciences tools provider with sales in over 100 countries. The group is active in lab equipment, pharma equipment and medical devices. Manufacturing facilities are located in Asia and Europe. R&D is conducted worldwide spanning the US, Europe and Asia. Sales, service and marketing subsidiaries are located in 12 major markets including the US, UK, Singapore, Japan, China and India. Regional distribution centers are located in the US, UK, and Singapore.

Life Science • Chemical Research • Assisted Reproductive
 Technology (ART) • Pharmaceutical Equipment • General Equipment

ESCO

WORLD CLASS. WORLDWIDE.

Esco Micro Pte. Ltd. • 21 Changi South Street 1 • Singapore 486 777
 Tel +65 6542 0833 • Fax +65 6542 6920 • mail@escoglobal.com
 www.escoglobal.com

Esco Technologies, Inc. • 903 Sheehy Drive, Suite F, Horsham, PA 19044, USA
 Tel +1 215-441-9661 • Fax 484-698-7757 • eti.admin@escoglobal.com
 www.escolifesciences.us

Esco Global Offices: Bangladesh | Cameroon | China | Denmark | Germany | Hong Kong | India | Indonesia | Italy | Japan | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | United Kingdom | USA | Vietnam

ISOCIDE™



9010960 Powder Weighing Balance FW1-A, Sell Sheet, A4, vD, 07/09/19
 Esco can accept no responsibility for possible errors in catalogues, brochures and other printed materials.
 Esco reserves the right to alter its products and specifications without notice. All trademarks and logos
 in this material are the property of Esco and the respective companies.